



IF IT'S WORTH BUILDING, IT'S WORTH TESTING

This wishlist was generated on 28/07/2017, and contains the following Products:

36-3090/01

ADR Touch 2000 Standard Compression Machine with Digital Readout.

37-6130

100kN Flexural (Beams) Frame. Supplied Without Specimen Bearers and Fitting Kit.

37-6138

100kN Flexural Fitting Kit (ADR) Used for Connecting Flexural Frames to ADR Compression Machines.

37-6131

BS 1881 Specimen Bearer Assembly 38mm Dia X 160mm Length.

37-6133

Ball Seating Assembly 150mm Dia Platens.



ADR Touch 2000 Standard Compression Machine with Digital Readout.

Code: [36-3090/01](#)

Product Group: [General Purpose Compression Machines](#)

- 2000 kN/450 000 lbf capacity
- Tests 150 and 100 mm concrete cubes or cylinders up to 320 x 160 mm diameter
- Supplied with Windows(R) download software as standard.

Incorporating the ADR Touch digital readout, the machines are designed to test cubes and cylinders in accordance with most International Standards. Supplied fitted for cylinder testing with safety gates. When used for cube testing appropriate distance pieces according to the size of specimen to be tested are required and must be ordered separately.

Further Information

****NOT AVAILABLE FOR SALE IN THE USA****

The Standard range of compression machines has been designed to meet the need for a simple, economic and reliable means of testing concrete.

Whilst delivering all the features and reputation of the established ADR Series with its extensive design history, the new and improved user interface provides a high quality platform for testing that enhances the performance of ELE's compression machines.

New, sophisticated electronics further the benefits of simplified operation, whilst delivering the highest levels of accuracy in testing concrete and cement/mortar samples, satisfying the needs of Quality Control Managers, Lab Managers and Technicians.

ACCURACY & SAVINGS

The new ADR Touch Series, with 145mm (5.7") high resolution QVGA touch screen interface and intuitive menu-driven operation, reduces the time taken to set up the machine and perform tests, reducing the time to train staff by up to 25%.

- Up to 6 sample types can be set as favorites, enabling one-touch set up for repeat testing.
- Full, QWERTY touch pad for input of test data.

Most testing errors produce lower strength results. Noncompliant loading rates can generate errors in measured strength. The user interface includes real-time display of load vs. time, further ensuring accurate and consistent test results and providing "goodness of test" data to improve traceability in your QC operations.

TRACEABILITY & DATA QUALITY

The ADR range now provides improved data quality and traceability in due diligence cases - it is now possible to demonstrate traceability all the way from the machine/user to the accreditation body, increasing your reputation and peace of mind - all test results now come complete with the machine serial number attached.

- Enhanced USB data communications between PC and machine - eliminating the need for download software.
- Two gigabytes of storage memory.
- Full customization of sample sizes - stress calculations are automatically recalculated.


USER SAFETY

With full safety gates as standard, total systems diagnostics, ram run-out switches and overload warnings ensure the safety of your employees and the reliability of your machine.

Specification

Testing Standards	BS EN ISO 7500-1, ASTM E4
Capacity	2000 kN
Cubes (Concrete)	Up to 150 mm
Cylinders (Concrete)	Up to 160 x 320 mm
Blocks	n/a
Flexural Testing	Via Flexural Frame
TFV and ACV	Yes
Voltage Supply	220-240V 50/60 Hz 1Ph.
Frame Type	Welded
Max vert. Clearance	340 mm
Max hor. Clearance	325 mm
Platen Sizes	Lower, Upper 222 mm
Max Ram Travel	50 mm

Accessories

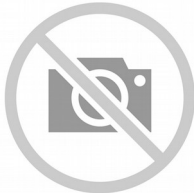

100kN Flexural Fitting Kit (ADR) Used for Connecting Flexural Frames to ADR Compression Machines.

Code: [37-6138](#)



2000 kN Capacity Calibration Load Cell with Hand-Held Readout

Code: [37-8315](#)



ADR Touch Head for Compression machines

Code: [37-4950/09](#)



Flexible Hose

Code: [5038B0027](#)



Impact Printer Rs232 Serial Connection. Supplied Complete with 1 Paper Roll.

Code: [37-4859/01](#)



Standard Distance Piece - 100 mm Depth

Code: [37-5100](#)



Standard Distance Piece - 20 mm Depth

Code: [37-4980](#)



Standard Distance Piece - 50 mm Depth

Code: [37-5000](#)



Standard Distance Piece - 60 mm Depth

Code: [37-5020](#)



Standard Distance Piece - 80 mm Depth

Code: [37-5050](#)

Spares/Consumables



Pressure transducer, 0-700bar, 0.05 to 10v DC Output

Code: [6014A0062](#)

Alternatives



ADR Touch 1500 Compression Machine with Digital Readout.

Code: [36-0720/01](#)



100kN Flexural (Beams) Frame. Supplied Without Specimen Bearers and Fitting Kit.

Code: [37-6130](#)

Product Group: [Flexural Testing for Beams](#)

- Open sided for ease of specimen loading
- Meets EN 12390-5, 1521, 13161, 772-6 and ASTM C78 for standard section beams
- Optional ball seating assembly

This rigidly constructed, open sided frame is suitable for testing 100 mm and 150 mm square-section beams. The frame supports a hydraulic ram and upper sub-platen assembly incorporating a spherical seating. The upper and lower sub-platens will accept various specimen loading bearers, which are supplied separately.

Standards

BS 1881-118, EN 12390-5

Further Information

Not available for sale in the USA

Accessories



100 kN Flexural Fitting Kit ADR Auto for 110 - 120 V AC, 60 Hz, 1 ph

Code: [37-6135/02](#)




100 kN Flexural Fitting Kit ADR Auto for 220 - 240 V AC, 50 Hz, 1 ph

Code: [37-6135/01](#)



100 kN Flexural Fitting Kit ADR Auto for 220 - 240 V AC, 60 Hz, 1 ph

Code: [37-6135/06](#)

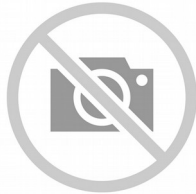

100kN Flexural Fitting Kit (ADR) Used for Connecting Flexural Frames to ADR Compression Machines.

Code: [37-6138](#)



ASTM C78 Specimen Bearer Assembly Case-Hardened 38mm Dia X 160mm Length.

Code: [37-6132](#)



Ball Seating Assembly 150mm Dia Platens.

Code: [37-6133](#)



BS 1881 Specimen Bearer Assembly 38mm Dia X 160mm Length.

Code: [37-6131](#)



Standard Distance Piece - 100 mm Depth

Code: [37-5100](#)



Standard Distance Piece - 20 mm Depth

Code: [37-4980](#)



Standard Distance Piece - 50 mm Depth

Code: [37-5000](#)



Standard Distance Piece - 60 mm Depth

Code: [37-5020](#)



Standard Distance Piece - 80 mm Depth

Code: [37-5050](#)

Alternatives



Accu-Tek Touch 250 Automatic Compression Machine 110 V 60Hz

Code: [36-1690/02](#)



Accu-Tek Touch 250 Automatic Compression Machine 220-240 V 50Hz

Code: [36-1690/01](#)



Accu-Tek Touch 250 Automatic Compression Machine 220-240 V 60Hz

Code: [36-1690/06](#)



ADR Touch Manual Power Pack Including Fitting Kit 110V 60Hz 1Ph

Code: [36-6165/02](#)



ADR Touch Manual Power Pack Including Fitting Kit 220-240V 50/60Hz 1Ph

Code: [36-6165/01](#)



100kN Flexural Fitting Kit (ADR) Used for Connecting Flexural Frames to ADR Compression Machines.

Code: [37-6138](#)

Product Group: [Accessories to connect ELE compression machines to the flexural/transverse frames](#)

These flexural frames can be used with a wide range of ELE compression machines. Simply select one of two fitting kits designed to link the flexural frame to either the ADR-Auto range of compression machines or any other ADR compression machine.

Further Information

Not available for sale in the USA



BS 1881 Specimen Bearer Assembly 38mm Dia X 160mm Length.

Code: [37-6131](#)

Product Group: [Flexural Testing for Beams](#)

Comprising 2 self-aligning upper bearers, 1 self-aligning and 1 fixed lower bearer. Roller bearers are 38 mm diameter x 160 mm long. Suitable for 3- or 4-point flexural testing of beams.

Standards

BS 1881-118, EN 12390-5

Further Information

****Not available for sale in the USA****



Ball Seating Assembly 150mm Dia Platens.

Code: [37-6133](#)

Product Group: [Flexural Testing for Beams](#)

Comprising a ball seating assembly with 150 mm diameter platen and a lower platen 150 mm diameter x 16 mm thick. Suitable for testing low-strength specimens.

Further Information

****Not available for sale in the USA****