

















IF IT'S WORTH BUILDING, IT'S WORTH TESTING

### **Contents**





Page 2 of 38

This wishlist was generated on 20/07/2018, and contains the following Products:

39-0045/01

**Automatic/Manual 5 litre nominal capacity Mortar Mixer** 

39-1100

**Three-gang Mould** 

24-9010

Straight Edge 300mm.

39-1150/01

**Jolting Table 220 - 240 V 50 Hz** 

39-1120

**Feeding Hopper** 

38-2015/01

**Automatic Vicat Apparatus** 

38-2015/12

**Needle Cleaning Device** 

38-2015/14

**Printer Paper Rolls** 

38-2015/15

**Mould Tank** 

38-2015/16

**Consistency Plunger** 

38-2021

**EN Initial and Final Set Needle** 

38-2300

**EN Vicat Mould** 

82-5270

Non Hazardous Glass Thermometers -1 - +51°C

38-3005

Le Chatelier Mould

www.ele.com # +44 (0)1525 249 200 # +1 (800) 323 1242

## **Contents**

www.ele.com +44 (0) 01525 249 200

38-3200

**Extensibility of Mould Apparatus BS** 

34-0140

300 mm Stainless Steel Rule.

38-3300

**Tamping Rod** 

38-3420/01

Le Chatelier Water Bath

39-1170

**Mortar Sand** 

82-3320

**Porcelain Crucible** 

39-7160/01

**ADR Touch Control PRO 250/25 Auto Cement Machine** 

39-6220

**Set of Distance Pieces** 







# Automatic/Manual 5 litre nominal capacity Mortar Mixer

Code: 39-0045/01

Product Group: Automatic/Manual 5 Litre Nominal Capacity Mortar Mixer

- New Microprocessor control
- New mix program selector
- Complies with latest EN Standards update
- Choice of automatic mixing cycles
- Sand and water dispenser supplied as standard

This mixer is designed to mix mortars and cement pastes to the requirements of the EN 196. The mixing paddle has a planetary motion and is driven by a motor with a microprocessor based speed and program controller. The mixer can be operated either in an automatic or manual mode. When the mixer is used in the manual mode, the two mixing speeds can be changed by means of a rocker switch, without switching off the motor. In the automatic mode any one of the preset mixing programs may be selected.

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

Complete with bowl and paddle, sand and Water dispensers. For 220V AC, 50 Hz, 1 ph.

#### **Specification**

Dimensions (l x w x h) 810 x 610 x 750 mm Speeds (rpm) Low, High Paddle:  $140 \pm 5$ ,  $285 \pm 10$ Mixing Head:  $62 \pm 5$ ,  $125 \pm 10$ 

Rated Power 180 W

Bowl Capacity 5 litres approx.

Weight 75 kg


www.ele.com +44 (0)1525 249 200

**=** +1 (800) 323 1242

Page 4 of 38





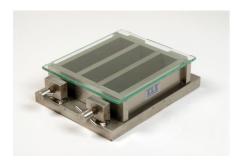

## **Three-gang Mould**

Code: 39-1100

Product Group: Moulding of Prisms

For  $40.1 \times 40 \times 160$  mm mortar prisms. This mould is manufactured to a very high specification. Supplied with glass plate.

#### **Standards**

EN 196-1, EN 413-2, EN 459-2, EN 1015-10, EN 13454-2, ISO 679, BS 3892-1, BS 4551-1

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

#### **Accessories**



#### **Feeding Hopper**

Code: 39-1120

### Spares/Consumables



#### **Glass Plate**

Code: 39-1100/10



#### **Mortar Sand**

Code: 39-1170

#### **Alternatives**

www.ele.com # +44 (0)1525 249 200 # +1 (800) 323 1242 Page 5 of 38

www.ele.com +44 (0) 01525 249 200





**Feeding Hopper** 

Code: 39-1120

www.ele.com +44 (0)1525 249 200







#### Straight Edge 300mm.

Code: 24-9010

Product Group: BS/EN compaction test 4.5kg, Proctor

compaction, Straightedge

300 mm

#### **Standards**

BS 1377-4, BS 1924-2, ASTM D698, ASTM D1557, AASHTO T99, EN 13286-2, EN 13286-4, ASTM D1883, AASHTO T193, ASTM D558, ASTM D559, ASTM D560, AASHTO T134, AASHTO T135, AASHTO T136, AASHTO T180

#### **Specification**

Construction Flat ground tool steel; beveled

edge.

Dimensions L x W x T (300mm x 40mm x

6mm).

Weight 550g

#### **Accessories**



#### **ASTM Compaction Rammer 2.5Kg.**

Code: 24-9063

#### **Alternatives**



#### **ASTM Compaction Rammer 2.5Kg.**

Code: 24-9063

www.ele.com # +44 (0)1525 249 200 # +1 (800) 323 1242 Page 7 of 38







### **Jolting Table 220 - 240 V 50 Hz**

Code: 39-1150/01

Product Group: Jolting Table

This machine consists of a mould table seated on a rotating cam driven at 60 revolutions per minute. The apparatus is supplied with a remote mains switch box, incorporating a push button start/stop control and automatic stop control at end of test.

Supplied without moulds.

#### **Standards**

EN 196-1, EN 413-2, EN 459-2, EN 1744-1, EN 1015-10, EN 1015-11, EN 13454-2, ISO 679, BS 3892-1, BS 4551-1

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

#### **Specification**

Motor rating 0.24 kW (1/3 hp)

Cam shaft RPM 60

Table lift 15 +/- 0.3 mm

Overall dimensions with swing (H x W x D)

bolts raised 340 x 400 x 1045 mm

Weight 90 kg

#### **Accessories**



#### **Feeding Hopper**

Code: 39-1120



#### **Three-gang Mould**

Code: 39-1100

www.ele.com # +44 (0)1525 249 200 # +1 (800) 323 1242 Page 8 of 38





## Spares/Consumables



#### **Mortar Sand**

Code: 39-1170



### Scraper

Code: 39-1130

#### **Alternatives**



Jolting Table 220 - 240 V AC 60 Hz

Code: 39-1150/06







## **Feeding Hopper**

Code: 39-1120

Product Group: Moulding of Prisms

For 39-1100 Mould.


#### **Standards**

EN 196-1, EN 413-2, EN 459-2, EN 1744-1, EN 1015-10, EN 13454-2, ISO 679, BS 3892-1

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

#### **Accessories**

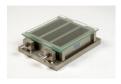

#### **Three-gang Mould**

Code: 39-1100

#### **Alternatives**



### **Three-gang Mould**

Code: 39-1100

www.ele.com +44 (0)1525 249 200







### **Automatic Vicat Apparatus**

Code: 38-2015/01

Product Group: Automatic Vicat

Complete with EN and ASTM Initial and Final needles, Consistency Plunger, Windows(R) software and RS232 cable, 1 x EN and ASTM Mould and Glass Plate. Suitable for continuous use in saturated humidity, at a controlled temperature of  $20^{\circ}\text{C} \pm 1^{\circ}\text{C}$ .

#### **Standards**

ASTM C187, ASTM C191


#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

For 220 - 240 V AC, 50 - 60 Hz, 1 ph.

#### **Accessories**

**ASTM Vicat Mould** 

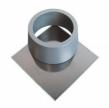

Code: 38-2660

www.ele.com # +44 (0)1525 249 200 # +1 (800) 323 1242 Page 11 of 38







**EN Vicat Mould** 

Code: 38-2300



#### **Mould Tank**

Code: 38-2015/15



## **Needle Cleaning Device**

Code: 38-2015/12



## **Printer Paper Rolls**

Code: 38-2015/14



#### **Vicat Mould**

Code: 38-2200

www.ele.com +44 (0)1525 249 200 = +1





## Spares/Consumables



#### **ASTM Initial and Final Set Needle**

Code: 38-2023



## **Consistency Plunger**

Code: 38-2015/16



#### **EN Initial and Final Set Needle**

Code: 38-2021

**= +1 (800) 323 1242** Page 13 of 38 www.ele.com







## **Needle Cleaning Device**

Code: 38-2015/12

Product Group: Automatic Vicat

Accessory.

#### **Standards**

ASTM C187, ASTM C191

### **Further Information**







### **Printer Paper Rolls**

Code: 38-2015/14

Product Group: Automatic Vicat

Accessory: pack of 10.

#### **Standards**

ASTM C187, ASTM C191, EN 196-3, EN 13454-2

#### **Further Information**







#### **Mould Tank**

Code: 38-2015/15

Product Group: Automatic Vicat

Accessory: for testing samples under water as per EN 196-3, for use in temperature controlled laboratories.

#### **Standards**

ASTM C187, ASTM C191

#### **Further Information**







## **Consistency Plunger**

Code: 38-2015/16

Product Group: Automatic Vicat

Spare.

#### **Standards**

ASTM C187, ASTM C191

#### **Further Information**





#### **EN Initial and Final Set Needle**

Code: 38-2021

Product Group: Automatic Vicat

#### **Standards**

ASTM C187, ASTM C191


#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

+44 (0)1525 249 200





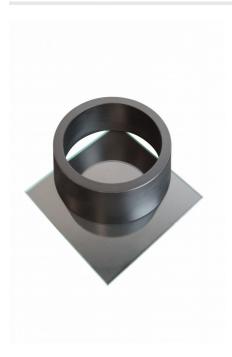

#### **EN Vicat Mould**

Code: 38-2300

Product Group: Vicat Method, Automatic Vicat

Accessory: manufactured from a hard rubber compound and supplied complete with a glass base

plate.

#### Standards

EN 196-3

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

#### **Accessories**


#### **ASTM Initial Set Needle**

Code: 38-2620







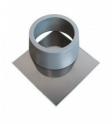

#### **ASTM Vicat Mould**

Code: 38-2660

#### **EN Initial and Final Set Needle**

Code: 38-2110

### **Vicat Mould**

Code: 38-2200



## Spares/Consumables

## **Consistency Plunger**

Code: 38-2010/10

#### **Alternatives**

www.ele.com # +44 (0)1525 249 200 # +1 (800) 323 1242 Page 20 of 38





Page 21 of 38


#### **ASTM Initial Set Needle**

Code: 38-2620



#### **ASTM Vicat Mould**

Code: 38-2660

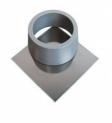

#### **EN Initial and Final Set Needle**

Code: 38-2110



#### **Vicat Mould**

Code: 38-2200



www.ele.com

**★** +44 (0)1525 249 200 **■** +1 (800) 323 1242







# Non Hazardous Glass Thermometers -1 - +51°C

Code: 82-5270

Product Group: Non Hazardous Glass Thermometers

0.1°C divisions + total immersion

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

www.ele.com +44 (0)1525 249 200







#### Le Chatelier Mould

Code: 38-3005

Product Group: General Expansion of Dry Hydrated Lime by the Le Chatelier Method, Soundness of Cement Paste by the Le Chatelier Method, Soundness

Comprising a split cylinder fitted with two indicator stems. The mould is supplied complete with two glass plates and a weight 100 g ±10 g. Three moulds required for each test.

#### **Standards**

BS 6463, EN 196-3, EN 459-2

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

#### **Accessories**



#### **Extensibility of Mould Apparatus BS**

Code: 38-3200



#### **Extensibility of Mould Apparatus EN**

Code: 38-3205



#### Le Chatelier Water Bath

Code: 38-3420/01

**= +1 (800) 323 1242** Page 23 of 38







**Tamping Rod** 

Code: 38-3300



**Trowel - Gauging** 

Code: 81-0335







## **Extensibility of Mould Apparatus BS**

Code: 38-3200

Product Group: General Expansion of Dry Hydrated Lime by the Le Chatelier Method

Conforms to BS 6463. The unit comprises a metal sleeve with a hook and set screw to fit over one of the mould pointers, and a clamp to fit onto the other pointer of the mould. Supplied complete with one weight 100 g  $\pm 1$  g.

#### **Standards**

EN 196-3

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

#### **Accessories**



#### Le Chatelier Water Bath

Code: 38-3420/01



#### **Tamping Rod**

Code: 38-3300



**Trowel - Gauging** 

Code: 81-0335

www.ele.com +44 (0)1525 249 200





#### **Alternatives**



**Le Chatelier Water Bath** 

Code: 38-3420/01



**Tamping Rod** 

Code: 38-3300



**Trowel - Gauging** 

Code: 81-0335







### 300 mm Stainless Steel Rule.

Code: 34-0140

Product Group: Slump Test

Standards EN 12350-2

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

www.ele.com +44 (0)1525 249 200







## **Tamping Rod**

Code: 38-3300

Product Group: General Expansion of Dry Hydrated

Lime by the Le Chatelier Method

17 mm end diameter.

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

#### **Accessories**



### **Extensibility of Mould Apparatus BS**

Code: 38-3200



### **Extensibility of Mould Apparatus EN**

Code: 38-3205



#### Le Chatelier Water Bath

Code: 38-3420/01

#### **Alternatives**

www.ele.com # +44 (0)1525 249 200 # +1 (800) 323 1242 Page 28 of 38





Page 29 of 38



### **Extensibility of Mould Apparatus BS**

Code: 38-3200



## **Extensibility of Mould Apparatus EN**

Code: 38-3205



**Le Chatelier Water Bath** 

Code: 38-3420/01







#### Le Chatelier Water Bath

Code: 38-3420/01

Product Group: Soundness of Cement Paste by the Le Chatelier Method, Soundness

Manufactured from corrosion resistant material, complete with a removable rack to hold up to twelve moulds and timer unit.

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

#### **Accessories**



#### **Extensibility of Mould Apparatus BS**

Code: 38-3200



### **Tamping Rod**

Code: 38-3300



**Trowel - Gauging** 

Code: 81-0335

#### **Alternatives**

www.ele.com # +44 (0)1525 249 200 # +1 (800) 323 1242 Page 30 of 38





Page 31 of 38



### **Extensibility of Mould Apparatus BS**

Code: 38-3200



## **Tamping Rod**

Code: 38-3300



**Trowel - Gauging** 

Code: 81-0335





#### **Mortar Sand**

Code: 39-1170

Product Group: Moulding of Prisms

Graded pack for 3 prisms. 1350g pack.

#### **Standards**

EN 196-1, EN 413-2, EN 459-2, EN 1744-1, EN 1015-10, EN 1015-11, EN 13454-2, ISO 679, BS 3892-1, BS 4551-1

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

#### **Accessories**



**Three-gang Mould** 

Code: 39-1100

www.ele.com +44 (0)1525 249 200







#### **Porcelain Crucible**

Code: 82-3320

Product Group: Loss-on-Ignition, Bottles and

Crucibles

30ml capacity complete with lid.

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

#### **Alternatives**



### **Heat of Hydration Apparatus**

Code: 38-4600/06



### **Heat of Hydration Apparatus**

Code: 38-4600/01







# ADR Touch Control PRO 250/25 Auto Cement Machine

Code: 39-7160/01

Product Group: Activity Index, ADR Touch Control Pro 250/25 Auto Cement Machines

- 250 kN maximum capacity
- Calibration accuracy to BS EN ISO 7500-1; ASTM
- Automatic loading cycle
- 25 kN low capacity frame supplied as standard
- Tests a wide variety of specimen sizes
- Tests mortar, lime, cement and Fly Ash
- Complete with compression/flexural jigs and platen sets
- Supplied with Windows(R) download software as standard

The ADR Touch Control PRO 250/25 Machine provides consistent automatic testing of a wide range of specimens. The machine comprises a standard ELE 250/25 kN load frame and the new ADR Touch Control PRO console and incorporates all the features included in the ADR Touch Control PRO series concrete testing machines.

As standard the machine is supplied with platens fitted to the load frame, compression jig with 40 mm and 50mm/2 inch square platen sets and flexural jig for testing  $40.1 \times 40 \times 160$  mm prisms.

The availability of the 25 kN low capacity load frame as standard extends the test capability of the machine for low strength compression or flexural testing.

#### **Standards**

EN 196-1, EN 459-2, EN 1744, EN 1015-11, EN 13454-2, BS 3892-1, BS 4551-1, ASTM C109, ISO 679

#### **Further Information**

\*\*Not available for sale in the USA\*\*

The ADR Touch Control PRO will deliver all the features and quality of the established ADR-Auto range, with its 20 year history, but with a new sleek design and additional capabilities. The console assembly consists of an ADR Touch Control PRO and power base which can be used with all existing concrete and cement frames

#### **USER INTERFACE**

The new ADR Touch Control PRO has an improved user interface providing a high quality platform

www.ele.com +44 (0) 01525 249 200



for testing that will enhance the performance of our compression machines. It now has a larger 7-inch colour touch screen with a higher definition display, allowing easier observation of test progress. A favourites option is also available, allowing users to save the most common set-ups for immediate access, without the need for navigating through selection menus.

#### **USER FLEXIBILITY**

The ADR Touch Control PRO has a unique flexible head with rotation capabilities of 75-degrees and a 45-degree tilt, improving usability and comfort for individual users. It is available with four languages English, Spanish, Portuguese and French, making it even more user-friendly and accessible to a wider audience. The new LAN connectivity feature lets the operator use the PC based APP (ELE Logger) to download test results away from the test area. User flexibility has also increased with the new ADR Touch Control Pro as it can now be operated remotely, so tests can be initiated away from the test area. Many issues can be resolved quickly by our service team through the new remote diagnostics feature and updating to the latest version of software has never been easier.

#### **COMPATIBILITY**

The all new dual sensor capabilities allow users to test both concrete and cement samples on one frame. The ADR Touch Control PRO is compatible with all existing compression frames and upgrade kits are also available. The upgrade kits can be used to convert existing ADR-Auto consoles to the new Automatic Digital Control Readout.

#### **TESTING TIMES**

The new generation ADR Compression Machines have an increased control over the pace rate, which allows users to run the pace rate higher than average, but still within standard. This produces an average time saving of 10% for a typical 150mm cube sample.

#### **ELECTRICAL**

Reduced component count and no wire design results in increased reliability and serviceability. New low voltage interlocks and e-stop has increased safety. The ADR Touch Control PRO has the latest Arm Core Processor technology which has improved the processing power and speed.

#### **Accessories**



40mm Square Platen Set

Code: 39-5600/10



Impact Printer Rs232 Serial Connection. Supplied Complete with 1 Paper Roll.

Code: 37-4859/01

www.ele.com +44 (0)1525 249 200





Page 36 of 38



### **Pair of Square Platens 2"**

Code: 39-5600/11



#### **Set of Distance Pieces**

Code: 39-6220

## Spares/Consumables



### **Compression Frame Jig Assembly**

Code: 39-5600



### **Flexural Jig Assembly**

Code: 39-6160

#### **Alternatives**



#### ADR Touch Control PRO 250/25 Auto Cement Machine

Code: 39-7160/06

www.ele.com +44 (0)1525 249 200 = +1 (800) 323 1242








# ADR Touch Manual with 250kN Frame and Jig 220-240v 50/60Hz

Code: 36-6166/01

www.ele.com 🚆 +44 (0)1525 249 200 📕 +1 (800) 323 1242 Page 37 of 38





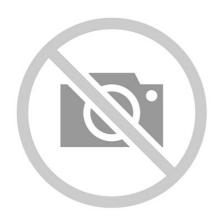

#### **Set of Distance Pieces**

Code: 39-6220

Product Group: Cement Testing Accessories

Accessory: for testing samples 40mm, 50mm, 70.7mm and 100mm. Note: These items are required when testing above sample sizes NOT using the compression jig.

#### **Further Information**

\*\*NOT AVAILABLE FOR SALE IN THE USA\*\*

+44 (0)1525 249 200

www.ele.com